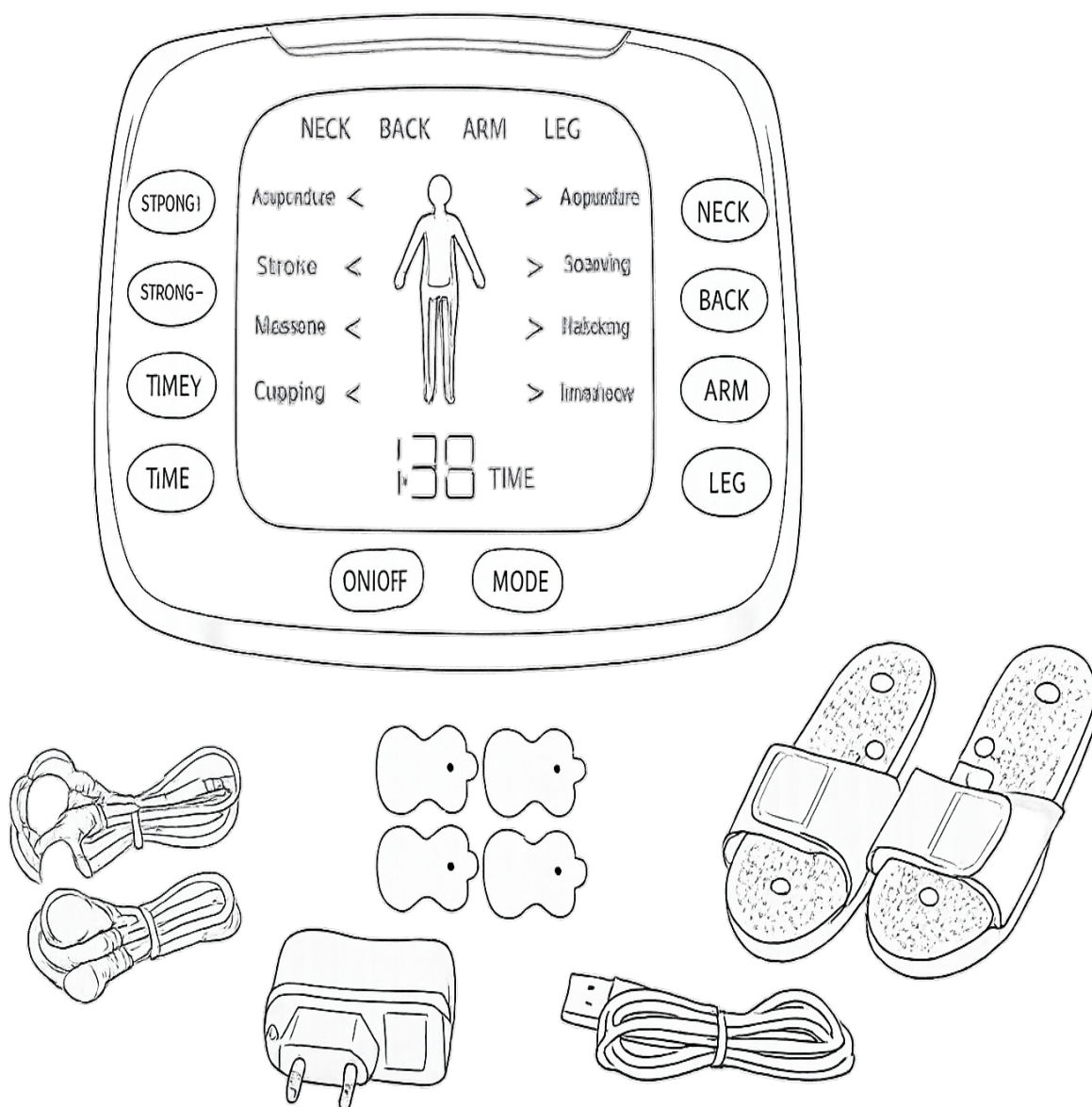
INSTRUKCJA

OBSŁUGI

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Dear Sir or Madam, thank you for purchasing our product!

Before using the product, please read carefully with the instructions below for proper use of the product.

Please keep this manual for future reference and following its recommendations, because failure to follow it may pose a threat to life or health.

Introduction to the user manual

This product is only suitable for decorative lighting of dry indoor rooms and is intended for private use.

This product is not intended for commercial or other uses. This product is not suitable for room lighting in the household.

Before first use, please read the instructions carefully. The information contained herein will also help you properly store the device, especially if it will not be used for an extended period.

Following the recommendations below will allow you to enjoy trouble-free and safe operation of the product.

ATTENTION!

Using the device in a manner inconsistent with the information contained in the manual may result in the risk of bodily injury.

Application and description of the device

The electromuscular stimulator is a device designed for home pulse therapy, supporting relaxation, regeneration, and muscle pain relief. It can be used on the neck, back, arms, and legs, as well as for muscle fatigue after physical exertion. The device offers several massage modes and adjustable intensity and duration, allowing you to tailor the experience to your individual needs.

This device is intended for home use only and is not intended for medical or commercial use, and is not intended to replace medical therapy. The device is intended for indoor use in dry rooms. It is not intended for medical, therapeutic, or commercial use.

For safety and CE certification reasons, this product must not be rebuilt or modified in any way. Using the product for purposes other than those described may result in damage to the product. Improper use may also result in hazards such as short circuits, fire, electric shock, etc.

WARNING!

- Danger of improper use!
- Using the device in a manner other than its intended purpose and/or using the device in any other way may involve various risks.

- The device should only be used for its intended purpose.
- Please follow the procedures described in this instruction manual.
- Keep the product away from children.
- Stop the massage if you experience discomfort, pain or dizziness.

TIP

- The device must be checked for completeness of delivery and any visible damage.
- In the event of an incomplete delivery or damage due to faulty packaging or transport, please contact the service hotline.
- It is best to do the massage for 10–15 minutes a day.

Technical data of the device

Symbol	24148
EAN	5907451335444
Material	plastic
Power supply	3x AA batteries or power supply
Flap insert length	26cm
Number of massage modes	8
Adjusting the massage length	Yes
Adjusting the massage intensity	Yes

Start-up / assembly instructions

- Note: Packaging material must be completely removed from the product.
- Insert the batteries as marked, paying attention to the correct polarity (+ and -) or connect the device to the power supply using the included power supply/cable.
- Connect the cables to the appropriate sockets on the device.
- Connect the electrodes to the wires
- Remove the protective foil from the electrodes and stick them to clean, dry skin in the desired location.
- Make sure the electrodes adhere evenly to the skin.
- Press the ON/OFF button to turn on the device.
- Select massage area (NECK, BACK, ARM, LEG).
- Select the operating mode using the MODE button.
- Set the pulse intensity using the STRONG+ / STRONG- buttons.
- Set the treatment time using the TIME+ / TIME- buttons
- After finishing the massage, turn off the device using the ON/OFF button.
- Carefully peel the electrodes from the skin and cover them with protective foil.

Selecting the electrical power supply

If the DC plug is not connected to the product, power is supplied by batteries.

Battery replacement instructions:

- Please insert the batteries correctly, paying attention to the +/- polarity.
- Only one type of battery should be used.
- Always replace the entire set of batteries.
- Exhausted batteries should be removed from the device and never thrown into fire.
- Never try to charge regular batteries.
- Store unused batteries in their packaging away from metal objects.
- **Do not leave children unsupervised.**

Safety tips and rules

- This product may be used by children aged 8 years and over and by persons with reduced physical, sensory or mental capabilities or lack of experience or knowledge, provided they are supervised or have been instructed on the safe use of the product and understand the resulting risks.
- Children are not allowed to play with the product.
- Cleaning and maintenance must not be performed by children without supervision.
- Never immerse the device in water.
- For cleaning, use a damp cloth or a mild detergent.
- Do not use a damaged device.
- Make sure packaging materials are not left unattended. Children may play with them, which is dangerous.
- Protect the product from extreme temperatures, direct sunlight, strong vibrations, high humidity, moisture, flammable gases, vapors and solvents.
- Do not expose the product to mechanical stress.
- If safe operation is no longer possible, discontinue use and secure the product against further use. Safe operation is no longer possible if the product: - has been damaged, - is not functioning properly, - has been stored for an extended period under unfavorable conditions, or - has been subjected to excessive stress during transport.
- Do not use the product if any part is damaged. If the cord is damaged, do not attempt to repair it yourself.
- Do not disassemble the device yourself.
- Protect the product from moisture.
- Lighting elements cannot be replaced. If the light elements stop working due to wear, the entire product must be replaced.
- Make sure not to damage the USB cable with sharp edges or hot objects. Fully unwind the USB cable before use.
- The included USB cable is only suitable for this product.
- Turn off the product immediately and remove the cable from the product if it smells burning or smoke is visible. Have the product checked by a qualified technician before using it again.
- If you're using a USB power adapter, the power outlet must be easily accessible at all times so that you can easily unplug the USB power adapter in an emergency. Please also follow the manufacturer's operating instructions.

- Do not use electrodes on damaged, irritated or infected skin.
- Do not place electrodes on the front of the neck or in the carotid sinus area.
- Do not use electrodes on the head or near the heart.
- Do not place electrodes on opposite sides of the chest.
- Do not use the device while bathing or showering.
- Do not use the device while sleeping.
- Do not drive or operate machinery during treatment.
- If you experience any discomfort, discontinue use immediately.
- After each use, clean your skin and cover the electrodes with protective foil.
- Do not use the device on people with a pacemaker or insulin pump.
- Do not use in pregnant or breastfeeding women.
- Do not use in people with cancer without consulting a doctor.
- If you have chronic diseases, consult your doctor before use.
- Do not exceed the recommended single treatment time (10 minutes).
- The product should always be used as intended.

Safety instructions for batteries/accumulators

- Keep batteries/rechargeable batteries out of the reach of children. If swallowed, seek medical attention immediately!
- Disposable batteries must not be recharged. Batteries/rechargeable batteries should not be short-circuited and/or opened. This may result in overheating, fire, or explosion.
- Never throw batteries/rechargeable batteries into fire or water.
- Batteries/rechargeable batteries should never be subjected to mechanical stress.
- Risk of acid leakage from batteries/accumulators.
- Avoid extreme conditions and temperatures that may affect batteries/accumulators, e.g. radiators/direct sunlight.
- If batteries/rechargeable batteries leak, avoid contact with skin, eyes, and mucous membranes! Immediately flush the affected areas with clean water and seek medical attention!
- Leaking or damaged batteries/rechargeable batteries can cause chemical burns when in contact with skin. Therefore, wear appropriate protective gloves.

Warning and information symbols

	<u>TIPS AND INFORMATION ON THE MANAGEMENT OF USED PACKAGING</u> The packaging is made of environmentally friendly materials that can be disposed of at your local recycling center. Used packaging material should be delivered to a designated waste disposal site designated by local authorities. For information on recycling options for used products, please contact your local municipality or city office.
---	---

	<u>DISPOSAL OF USED ELECTRICAL AND ELECTRONIC EQUIPMENT</u> For environmental reasons, used electrical and electronic products should not be disposed of with regular household waste, but should be disposed of properly. Information on collection points and their opening hours can be obtained from the relevant authorities.
	This product complies with the requirements of the applicable European and national directives. The product meets European and national requirements for the safety of devices and products.
	This product complies with the requirements of the relevant European and national RoHS directives.
	Improper disposal of batteries/accumulators poses a threat to the environment! Batteries/accumulators should not be disposed of with household waste. They may contain harmful heavy metals and should be treated as special waste. The chemical symbols for heavy metals are: Cd = cadmium, Hg = mercury, Pb = lead. Therefore, used batteries/accumulators should be taken to municipal hazardous waste collection points. The built-in battery must not be disassembled for disposal. The product should be taken to a collection point for used electronics.

We reserve the right to make changes to the text, design and product data without notice.

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